



Procedure and Examination Chairs

ACCESSORY CATALOG



Models

| S2 | S3 |

Sonesta manufactures and sells tables and chairs for urology, gynecology, urodynamic diagnostics and upper GI procedures.

Accessories: Next Generation Sonesta S2 & S3



Leg & foot rests

S2 S3 Art. number

Leg rest holder (single)

This leg rest holder is standard on all models. By loosening one handle, you can adjust the leg height, width or angle. Can be rotated 360 degrees.



S S

525-6338

Leg rests (pair)

This leg rest has a ball and socket construction allowing tilt of the leg cushion in any direction. In addition, it has infinite angling possibilities when mounted in our leg rest holder model 6338. Made of stainless steel.



S S

525-6336

Fold up footrests (pair)

These footrests are adjustable in length from 30–42 cm/11,8”–16,5”. They are also adjustable in width and slide outward 19 cm/7,5” each. The bottom is in the same upholstery as your chair. When folded they function as an extension cushion that transforms your chair into a table where the patient can lie down. The gas piston easily adjusts the angle. They are easily removed with a click-off system.



A S

525-6331

Footplates (single)

Made of stainless steel with a plastic foot section, fits into our standard leg rest holder.



A A

524-6337

Stirrups and clamps (pair)

Stirrups that provide safe and quick positioning which improves surgical site access. The stirrups comes with silicone straps for easier cleaning and attachment. Lithotomy range: +85° to +55°, abduction range: +25° to -9°, boot rotational angle: +/- 25°.



A A

Stirrups and clamps 525-3100

Patient weight capacity: 350 lbs / 159 kg

Benefits & use:

- Multiple positioning options
- Patient comfort
- Enables optimal examination access
- Engineered ergonomics

S

Standard configuraton

A

Optional accessory

Accessories: Next Generation Sonesta S2 & S3



Seat cushions



Art. number

Seat Cushion with 170 mm cutout

A 'quick change' accessory seat cushion, with a 170 mm/6,7" deep cutout, that can be switched in seconds with any other seat cushion. You simply lift the front edge, remove it and replace with another type of cushion.



525-6366-3

Seat cushion with 220 mm cutout

A 'quick change' accessory seat cushion, with a 220 mm/8,7" deep cutout, that can be switched in seconds with any other seat cushion. You simply lift the front edge, remove it and replace with another type of cushion.



525-6366-4

Seat cushion with 80 mm cutout

A 'quick change' accessory seat cushion for pediatric care, with a 80 mm/3,1" deep cutout, that can be switched in seconds with any other seat cushion. You simply lift the front edge, remove it and replace with another type of cushion.



525-6366-5

Seat cushion without cutout

A 'quick change' accessory seat cushion that can be switched in seconds with any of our other seat cushions. You simply lift the front edge, remove it and replace with another type of cushion.



525-6365

Benefits & use:

- Multiple seating options
- Patient comfort
- Enables optimal examination access
- Easy to use and interchangeable
- Ergonomic and low profile cushions



Standard configuration



Optional accessory

Accessories: Next Generation Sonesta S2 & S3



Extension & neck cushions



Art. number

Extension cushion

"Slide-in" extension cushion fits into the two square pipes under the front of the seat cushion and is easily removed when the patient's legs are put onto the leg rests. Makes the chair into a 193 cm/76" long flat examination table where your patient can lie down comfortably.



525-6364

Puzzle piece extension cushion

Puzzle piece "slide-in" extension cushion that fills the standard 170 mm cushion cutout. Fits into the two square pipes under the front of the seat cushion and is easily removed when patient's legs are put onto the leg rests. Makes the chair into a 193 cm/76" long flat examination table where your patient can lie down comfortably..



525-6362

Neck pillow

A neck pillow, adjustable in length, in the same upholstery as your table. A strip attached to the paper holder allows you to quickly change its position for your patient's height.



525-6346

Benefits & use:

- Multiple positioning options
- Ability to transform your chair into a flat bed
- Patient comfort
- Adjustable in length, 143 cm—193 cm/56"—76"
- User-friendly
- Ergonomic and low profile cushions



Standard configuraton



Optional accessory

Accessories: Next Generation Sonesta S2 & S3



Funnels



Art. number

Funnel kit

A removable plastic funnel (top diameter 236 mm/9") with drainage. It is mounted on a fixed height 'quick change' holder.



524-6325-5

Drainage receptacle

Standard receptacle 60 mm/2,4" deep – A stainless steel receptacle with an adjustable 'quick change' holder.



524-6325-1

Funnel and cysto pan

The Cysto pan comes standard with a tray that mounts to the chair that holds the pan in place. The pan comes equipped with a screen, a plug and a hose that extends to 193 cm/76".



524-6325-10

Benefits & use:

- Practical use
- Functional design
- Easy to use and interchangeable



Standard configurator



Optional accessory

Accessories: Next Generation Sonesta S2 & S3



Controls

S2 S3 Art. number

Two motor foot control

Waterproof foot control for two motor chairs.



S -

524-6353

Three motor foot control

Waterproof foot control for three motor chair.W



- S

524-6357

Two motor hand control

Waterproof hand control for two motor chairs.



A -

524-P6411

Three motor hand control

Waterproof hand control for three mWotor chairs.



- A

524-6355

Benefits & use:

- Engineered ergonomics
- User-friendly
- Precise positioning

S

Standard configuraton

A

Optional accessory

Accessories: Next Generation Sonesta S2 & S3



Other accessories



Art. number

IV pole

An IV pole with holder that attaches to a support under the seat cushion. The IV pole will stay in the height relative to the patient as it moves as the table is adjusted. To be mounted on patient's right side of the chair.



524-6316

Arm rest (single)

The Sonesta arm rest is designed to be both an arm support and a handgrip, giving the patient many possible places to hold onto. Easily unfolded to allow transfer from a wheelchair.



525-6339-1 (L)
525-6339-2 (R)

Paper towel holder


The convenient placement of the paper towel holder enables immediate access. The paper towel holder is also combined with the manoeuvring handle for ultimate function.



525-6326

Restraint strap

Velcro straps for patient fixation. Pair of 100 cm/39" long and 10 cm/3,9" wide straps.

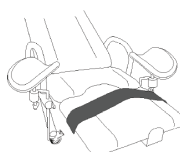

524-P6415

Benefits & use:

- Versatility
- Helpful examination tools
- Easy access



Standard configuraton



Optional accessory

How to accessorize after your and your patient's needs



Arm rests

Easy access to the patient from each side.



Leg rests

Swivelling leg rests for ultimate patient positioning and comfort.



Foot rests

Hydraulic piston helps switch from foot rests to extension cushion.



Foot plates

Convenient foot plates for function and comfort.



Neck pillow

An ergonomic accessory that enhances the patient comfort.



Puzzle piece

An extension cushion that transforms your chair into a table or flat bed.



Seat cushion

Choose the perfect fit for your patient and examination.



Motor control

Ergonomic motor controls, enables precise positioning.



Funnel

A functional funnel kit for performing non-UDS urology procedures.



Paper roll holder

Convenient placement for immediate access the paper towel roll.



Drainage receptacle

A stainless steel receptacle. Perfect for your instruments.



IV pole

Allows the patient to remain in one place during procedures.

An ADDvise Group Company

SONESTA

SONESTA MEDICAL AB
INDUSTRIVÄGEN 7
171 48 SOLNA, SWEDEN
+46 8 502 572 80
INFO@SONESTAMEDICAL.SE
ORDER@SONESTAMEDICAL.SE
WWW.SONESTAMEDICAL.SE

SONESTA MEDICAL INC.
2 DEBUSH AVENUE UNIT C3
MIDDLETON, MA 01949, USA
+1 630-519-3450
INFOUS@SONESTAMEDICAL.SE
ORDERUS@SONESTAMEDICAL.SE
WWW.SONESTAMEDICAL.SE